

## UNIT DOSE MEDICATION - SIMPLE TROLLEY 100 SERIES 1808 W - 10 PATIENTS

**DIMENSIONS** 640 x 480 mm

**ADEQUATE HEIGHT, ERGONOMICS** 1005 mm (work surface – floor).

### WORK SURFACE

- Polyethylene, PE. Asepsis, hygiene, easy-to-clean.
- Withstands alcohol, detergents, disinfectants, chlorhexidine.
- Protection on all 4 sides. Pleasant to the touch, rounded shapes.

### MONOBLOC STRUCTURE

- Stainless steel AISI-304. Fine satin finish, asepsis, hygiene.

### ISO DRAWER FRAMEWORK


- Total opening.
- Steel ball-bearing slides with self-closing system.
- Front stop, it protects the slide against strong impacts when closing.
- Load capacity per drawer: 20 kg.
- Interior frame ISO 600 x 400 mm.
- Polyethylene drawer front part, thickness 10 mm, 6 colours (Ref. G, Y, W, B, F, R).
- Ergonomic handle.
- Card holder with white cards, removable. ISO norm 7810.

### ISO DRAWERS 600 x 400 mm

- Heights 50, 100, 200 mm. Interchangeable, ABS.
- ISO norm 3394.

### CASSETTES FOR MEDICATION, INTERCHANGEABLE, 5 PATIENTS

- Width: 600 mm. ISO pharmacy furniture-compatible.
- Cassette resistant to daily use, made of polypropylene and epoxy painted steel, grey colour.
- Cassette locked. Manual tilting. Epoxy painted steel, grey colour.  
Prevents the access to medication to noauthorized people.  
Avoids the involuntary opening of bins when the cassettes are replaced from the transport trolley 600 series to a hospitalization units trolley 100 series.
- Medication bins, one for each patient. Grey polypropylene.  
Wide bin, designed for organizational needs of pharmacy in general hospitals.

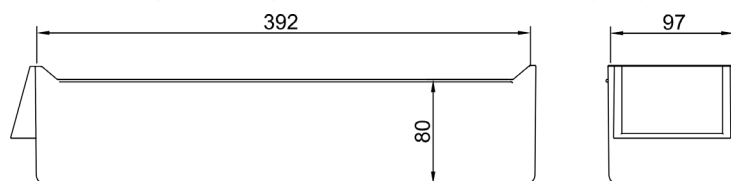

- Patient identity label, white, removable. ISO norm 7810.  
Cassette locked, upper part of the label hidden, lower part visible.
- Optional dividers, transparent polypropylene.

### KEY LOCKING

- Individual or multiple lock.
- Total locking of ISO drawers and cassettes for medication.
- Prevents from opening ISO drawers, the medication bins and the cassettes extraction.

### EXTENSIVE MOBILITY TROLLEYS, DAILY CONDUCTION

- Elevated push handle. Ergonomics, height 1045 mm.
- Double wheels Ø 125 mm. Soft tread, silent.
- 2 antistatic (electrical conductive), 2 with brake at the front part.

**1808 W - 10 PATIENTS**  
**108M6 W.32R.46.50**

**ACCESSORIES**

**32R** SIDE SUPPORT

**46** 3 GLOVES BOX HOLDER

**50** 13 LITRES BIN

